
Your reliable partner for APIs for over 70 years





Reliable supply chain



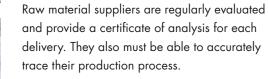





00% on-line inspection.

Guaranteed cleanliness and surface treatment

During washing and anodisation, the quality of ingredients and parameters is controlled by computer. The handling of containers is also fully automated. During several years, evaluation of the efficiency of the washing process has been performed and a log 3 (1/1000) efficiency factor has been demonstrated.

Packing and palletisation

Fully automated line from mechanical transformation until washing and anodisation. 100% control of the containers in a clean room with laminar flow before packing.

Batch and Unit traceability

Batch and unit traceability are ensured through bar codes labels and ink jet printing or engraving on the containers.

TOURNAIRE: Certified quality

Tournaire products meet today's most demanding quality standards and are manufactured in accordance with the strictest specifications.



EQUIPMENT I PACKAGING

WWW.TOURNAIRE.FR







FOURNAIRE Produces the best, protects the most

TOURNAIRE and APIs:

The ultimate protection for your active ingredients





Tournaire has registered since 2000 a Type III-C DMF with the US FDA which is updated and enriched annually.







TOURNAIRE Packaging: Major Partner for API containers since the 1950s

From streptomycin to last generations of cephalosporins, Tournaire has always provided state of the art containers for sterile or non sterile APIs. Our family owned company closely supports its clients with anticipation of new regulations (DMF), and mastership of modern techniques in the frame of the most demanding quality referentials.

Key figures:

- Over **180** years of experience.
- 300 staff in France.
- Over 40 000 m² of industrial buildings.
- 5000 tons of materials transformed
- Over **60** million € sales.
- 65% Export worldwide through an international network of agents and distributors.
- 2 manufacturing sites.

Tournaire Equipment:



Tournaire Equipment Division designs and manufactures production units for pharmaceuticals, flavors & fragrances in the field of extraction, evaporation, rectification, short path distillation, filtration, drying and distillation.

Its engineering department submits the technical files for their working and validation. Tournaire Equipment provides from pilot units to industrial plants compliant with ATEX and DESP rules.



Aluminium for stability and protection

Easy Sterilisation Depyrogenation

Dry heat sterilisation allows in one single step the sterilisation AND depyrogenation of the container. Furthermore, it eliminates the presence of moisture. Autoclaving can also be used for sterilisation alone. When tested before despatch, our containers have been found to be pyrogen free.



High resistance

Seamless aluminium containers are resistant to shocks, perforations, cuts, cold and heat. They also give absolute protection to light and gas permeation.

Safe and reliable protection

The European Pharmacopeia considers pure (min 99.5%) aluminium containers as safe and neutral for contact with Food and Drugs when produced and used according to GMP (Good manufacturing Practice).

Specific Surface treatment

Scouring and anodisation according to an approved process give a clean, smooth and resistant surface. Tournaire has worked closely with major pharmaceutical users over the years to standardise the best surface treatment, described in our DMF.



Washing

TOURNAIRE bottles:

Full range with type III-C DMF to help your registration



VOL.(L)	T802	T835	O 220
2,5	V	V	V
5,5	V	V	-
6,25	V	-	V
12,5	V	-	v
20	V	V	-
25	V	-	V
28,5	V	-	-
29,5	-	V	-
30	-	-	V
32	V	-	-

Geometrical features



Efficient design for user friendly solution

Seamless container

Each container is made of one disc of pure aluminium. Its whole surface is continuous, smooth and has no recesses nor

Optimised tightness

The calibration of the neck together with the elasticity of the rubber plug provides an hermetic seal, even with extreme variation of pressure.



Filling

Easy sealing and opening

Different closure systems are available, which provide user friendly sealing and opening as well as efficient tamper evident



TOURNAIRE WORKS FOR YOU by providing close regulatory

support, but also accessories in addition to its container range. Various validated closing systems and accessories will bring smooth and reliable use for your comfort and safety.

Easy to use closing systems

102mm opening with chlorobutyl plug and aluminium cap for sealing. 135 mm opening with chlorobutyl plug and aluminium cap for sealing.



220mm opening with aluminium lid, nitrile ring wad and clamping ring.



Sealing systems



steel crimping tool provides reinforced crimping efficiency sterilise and does not require maintenance nor adjustment. An automation of this manual

Cap Opener



Stainless steel cap opener allows to smoothly open the aluminium cap without touching it, thus avoiding the risk of contamination or cutting the glove.

Tool for inserting/ removing the plug



Connected with reduced vacuum source, the stainless steel removal/insertion tool allows to insert/remove the chlorobutyl plug easily (no rubbing nor effort) and avoids any contact with the plug (reduces risk of contamination).

Sample bottles



for 10 g to 100 g sampling.

